
FeelTrain1

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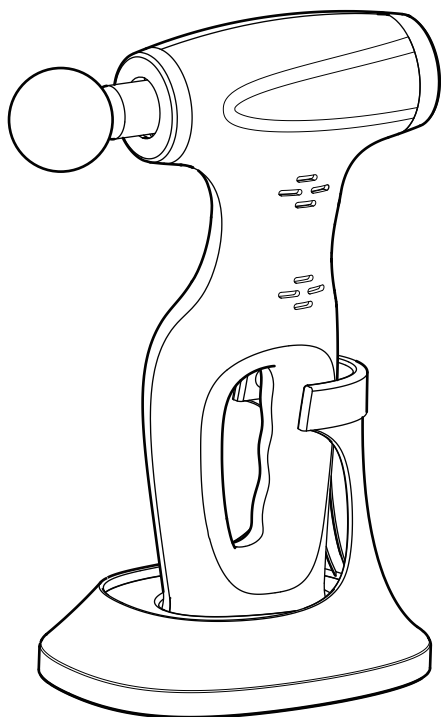Адрес производителя: D Building No. 100, JINFU ROAD, Chengnan Industrial Zone, TONG'AN District, Xiamen, China

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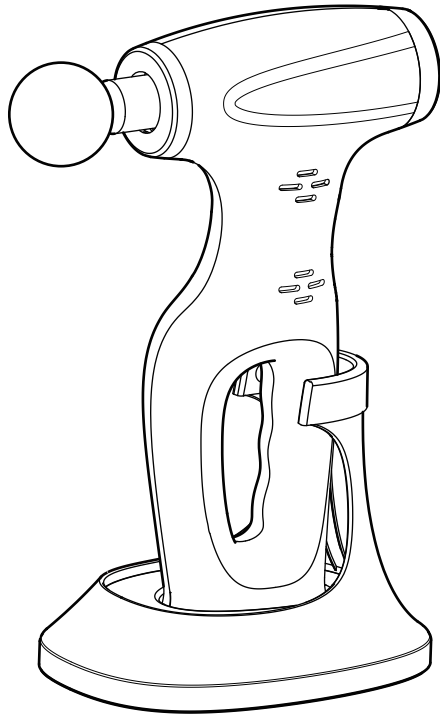
User Manual

Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER - To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Never operate the product with the air openings blocked. Keep the air openings away from lint, hair, and something like that.
- 4) Do not drop or insert any object into any opening.
- 5) Do not use outdoors.
- 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
- 8) Do not carry this product by its supply cord or use cord as a handle.
- 9) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
- 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
- 11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examinatio.
- 12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
- 13) Use this appliance only for its intended use as described in this manual.
- 14) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 15) The appliance is only to be used with the power supply unit provided with the appliance.
- 16) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

※ Children should be supervised to ensure that they do not play with the appliance.

※ Cleaning and user maintenance shall not be made by children without supervision.

17) This appliance contains batteries that are non-replaceable.

18) This appliance contains batteries that are only replaceable by skilled persons.

19) For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

* If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.

2. People who are supposed to lose their sense of temperature.

3. People with malignant tumors.

4. People with heart disease.

5. A person who is restless early in pregnancy or shortly after giving birth.

6. One who has an abnormal or curved spine.

7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.

8. Someone with a scar on their skin.

9. People who need quiet and people who are in poor health.

10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.

11. People with 38 °C temperature (fever).

12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.

13. People with tendonitis.

14. Others treated at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

• Do not wash!

• Do not dry clean!

Introduction

Thank you for purchasing the Gun Messenger. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and order number.

Product Specifications

Model	FeelTrain1 (Series JA)
Name	Gun Messenger
Input	DC17.2V === 1A
Default timer	10±1 minutes
Battery Specifications	DC 14.8V === 2200mAh


This product is intended for house-hold use only.

Short time operation 10 minutes, resting time 30 minutes

Note: Certified manufacturer's adapter must be equipped in the use of this products.

Package Includes

- 1 x Gun Messenger
- 5 x Massage Heads
- 1 x Adapter
- 1 x User manual
- 1 x Charging Dock

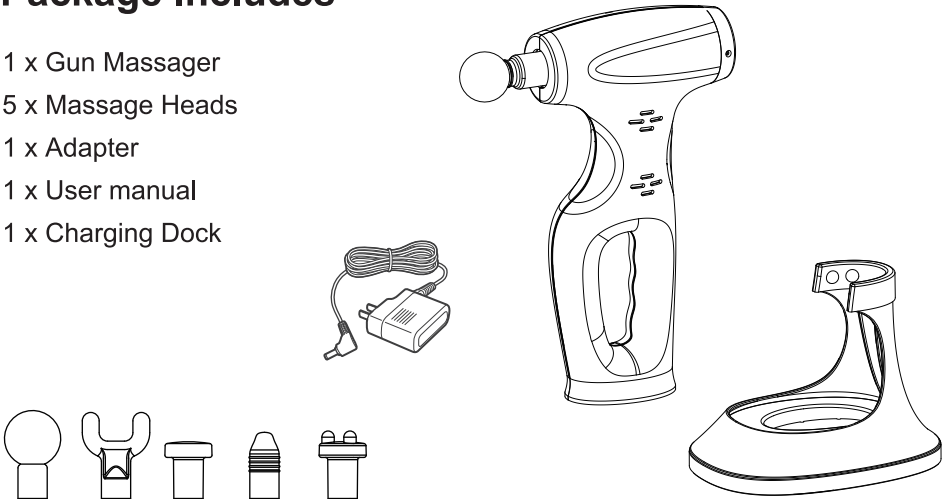
Product Operation Guide

Powering on Massager

- Select the proper massage head according to the actual situation and insert the massage head into the machine by rotating slightly.
- To power on your massager simply long press the  button for 3 seconds. Once that has been done, massager enters the 1st level intensity, 1 LED light on.
- Press the button again, massager switches to 2nd level, 2 LED lights on.
- Press the button again, turn on 3rd level, 3 LED lights on.
- Press the power button for the fourth time to back to the 1st level intensity, 1 LED light on...Loop in turn.

Powering off Massager


- In the process of starting up, long press the button for 3 seconds to turn off the massager at any intensity level.

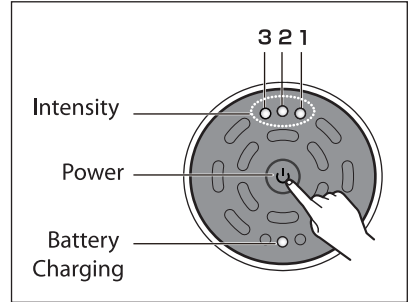Charging the Battery

- ① Place the massager into the charging dock.
 - ② Plug the power cable into the power port.
 - ③ Connect the adapter to a power source
 - ④ The power indicator lights up in red color when charging has started.
- The power indicator turns into green automatically when fully charged after 2 hours.

Note: When the battery is low and the massager is turned on, the power indicator will flash slowly in red color. Do not turn on the massager when it's being charged in the dock.

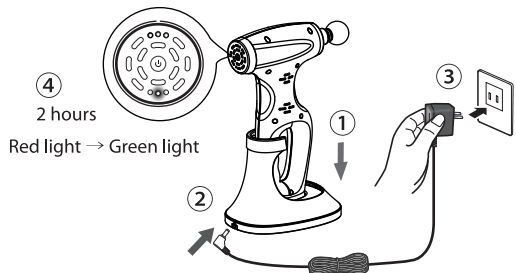
Removing the built - in rechargeable battery

Remove the built - in battery when disposing the appliance.
Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing the appliance, and must not be used to repair it.

If you dismantle the appliance by yourself, it may cause malfunction.

- Unplug the AC adaptor from the appliance before disassembling it.
- Turn on the power and then keep the power on until the battery is completely discharged.
- Perform steps ① to ② and lift the battery, and then remove it.
- Please take care not to short-circuit the battery.
- the appliance is only to be used with the power supply unit provided with the appliance.

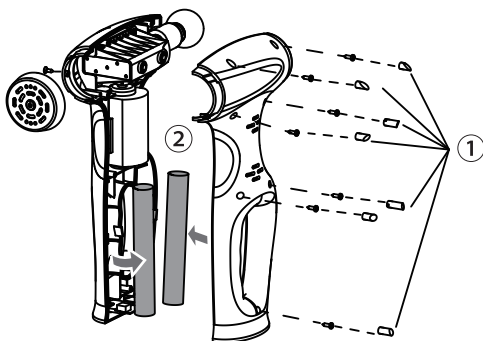How to clean

Cleaning

Before you clean, please make sure the power is switched off, and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.


Maintenance

Keep the product clean and far from the wet environment.
Clean when you feel cleaning is necessary.

Cleaning the main unit

- Wipe away dirt with a cloth soaked in soapy water, wring out the cloth before wiping.

Do not use alcohol, paint thinner, benzine, or nail polish remover. Using these may damage the appliance or cause cracks or discoloration in the components.

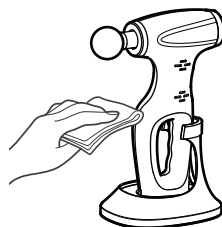Troubleshooting

Problem	Possible cause	Action
The appliance does not operate even the "power" switch is pressed.	▶ It is not charged	▶ Recharge with the AC adaptor. (See page 4.)
	▶ The "power" switch is not firmly pressed.	▶ Press the "power" button firmly to its place.
The charging indicator lamp does not go on/goes out.	▶ If the indicator lamp always go on when appliance is charging, that does not indicates a problem (See page 4)	
	▶ The adaptor is not properly inserted into the household outlet.	▶ Properly insert the adaptor into the household outlet.
	▶ The appliance plug is not firmly inserted into the appliance socket.	▶ Firmly insert the appliance plug into the appliance socket.
The appliance can be used only for few minutes even after charging it.	▶ Battery life is exhausted.	▶ Consult a retailer or an authorized service center. (Do not change the battery by yourself. Doing so may damage the appliance.)
The appliance stops during use.	▶ The battery has run down.	▶ Recharge with the AC adaptor.
	▶ The appliance automatically stops after approximately 10 minutes. This does not indicate a problem.	

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorised by us for repair.

	<h2>Correct Disposal of this product</h2>
	<p>This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.</p>

Technological parameter

Product: Gun Massager

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Input: DC17.2V === 1A

Battery Specifications: DC 14.8V === 2200mAh

Manufacturer: XIAMEN EMOKA HEALTHSCIENCE &
TECHNOLOGY CO. LTD.

Production date: on the massager

Made in China

Manufacturer's address: D Building No. 100, JINFU ROAD,
Chengnan Industrial Zone, TONG'AN District, Xiamen, China

Phone: 0086-755-25122982

General importer in Russia:

OOO "Max Concept", OGRN 1147453008532

Moscow, Lyubertsy district, Kotelniki, 1st Pokrovsky

PR-d, 5

phone 8-800-700-14-29

www.ergonova.ru

info@ergonova.ru

Full list of representative offices: <http://ergonova.ru/about>

/stores/ all copyrights to texts and images

owned by Ergonova GmbH and may not be used otherwise

only with the consent of the copyright holder.

Copyright © 2012.

Ergonova GmbH

All rights reserved.